humantouch.com
800.355.2762

©2010 Human Touch®, LLC. Patents pending and intellectual property rights enforced.

No medical claims are warranted or implied by the use of this product.


Human Touch is a trademark of Human Touch®, LLC.

HT-1280-NA-B0

ht 1280

use & care manual

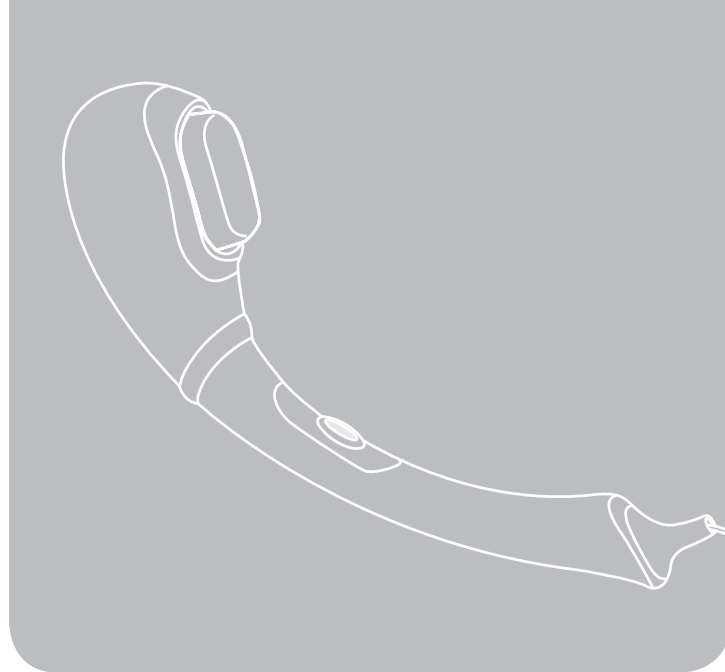ENG

use & care manual

FRC

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ESN

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ht 1280

welcome to **human touch**®

Your Human Touch® personal massager makes it easy to feel better whenever and wherever you'd like – its compact, design makes it perfectly portable. Your massager's unique shape makes it easy to target neck, shoulder, lower back, and other hard-to-reach muscles with deep-muscle relief. Loosen stiff muscles in the morning, or ease away stress and tension at the end of a long day. However you use it, your massager is a perfect way to rejuvenate.

Before using your massager, be sure to read this manual carefully.

Human Touch® **Partners in Your Well-Being**

For decades, Human Touch has been a leader in the field of massage products. We create state-of-the-art products that embrace the latest technology while applying the ancient healing art of massage. Balancing your physical and mental needs requires the right mix of proper nutrition, exercise, and massage. Human Touch® massage systems are vital components for achieving balance between physical and mental well-being.

To learn more about our company and our products, we invite you to visit us on the Internet at: www.humantouch.com or call our Customer Service Representatives at 800-355-2762.

Please send any questions or comments regarding this manual to documentation@humantouch.com

keep it safe

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before using this appliance.

 **DANGER: To reduce the risk of electrical shock:**

- Always unplug this product from the electrical outlet immediately after using and before cleaning.
- Do not reach for this product if it has fallen into water. Unplug immediately.
- Do not use while bathing or in the shower.
- Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place into or drop into water or other liquid.

 **WARNING: To reduce the risk of burns, fire, electric shock or injury to persons:**

- Do not use on sensitive skin areas or in the presence of poor circulation.
- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.

- Close supervision is necessary when this product is used by, on, or near children, invalids, or disabled persons.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return this product to a service center for examination and repair.
- Do not carry this appliance by supply cord or use cord as a handle.
- Always keep the power cord away from heated surfaces.
- Never operate this product with the air openings blocked. Keep the air openings free of lint, hair, and the like. Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
- Never drop or insert any object into any opening.
- While using the massager, keep your hair away from the air inlets.
- Do not leave unattended on any surface while it is operating.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- To disconnect, turn all controls to the OFF position, then remove plug from outlet.
- Do not massage abdomen when pain is present or if you are pregnant.
- Do not massage any area of the body that is swollen, inflamed or covered with skin eruptions.
- Do not use on unexplained pain.
- Do not use this product when standing on a damp floor, or when any part of the body is in contact with plumbing or any similar ground.

- Do not use this product if noise above normal sound is heard.
- Never use while sleeping.
- Always avoid contact with the eye or eyelid.
- Never recharge unit using any charger other than the charger that came with the unit.
- Never expose batteries to excessive heat.

Safety Precautions

- Comfortable clothing is recommended while using this product.
- If you feel any abnormal discomfort during the operation of this product, turn off the power immediately and do not use the product. Consult your physician.
- When using this product after it has been in storage, check to see that the product is operating properly prior to use.

SAVE THESE INSTRUCTIONS

Do not attempt to repair this product yourself.

This product is designed for personal in-home use. Commercial use will void warranty.

WARNING: To reduce the risk of electric shock, this appliance has a polarized plug (one blade is wide than the other). This plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. Do not change the plug in any way.

Care and Maintenance


If your massager becomes soiled, unplug the unit and wipe it clean with a cloth dampened with a non-abrasive cleanser. Do not use thinner, benzene or any other solvents. The product is specially designed and constructed to be maintenance free. No periodic lubrication or servicing is required. Your product should only be serviced by an authorized service center. Care should be taken not to overload the unit.

STORAGE:

Your product has been carefully designed to give you many years of trouble-free service. Always unplug the unit when not in use. Additionally, if not in use for a long period of time, we recommend that the power cord be coiled and that the unit be placed in an environment free of dust and moisture. Do not store near heat or open flame. Do not leave this product exposed to direct sunlight for long periods of time. This may cause fading or damage.

take control

1. Plug the massager into a power source.
2. Slide the power switch to the Slow (I) or Fast (II) speed position, depending on the desired massage intensity.
3. Use the massage head to soothe muscle stiffness away. For extra massage benefit, gently stretch the muscle as you are massaging it. It's best to use the massager on bare skin.
CAUTION: Do not massage any single area for more than 15 minutes. Do not massage the upper neck area, near the base of the skull.
4. When you have completed your massage. Slide the power switch to the off position (halfway between the I and II positions), then unplug the massager.

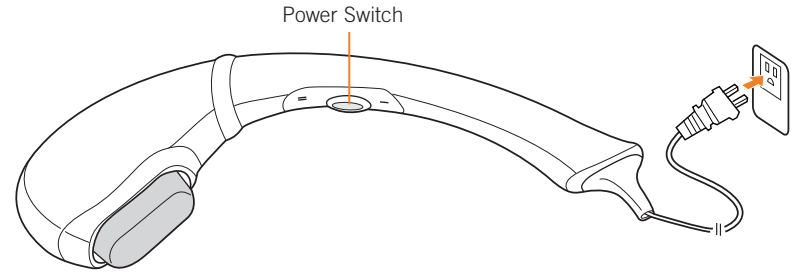Specifications*

Voltage: AC110-120V 60Hz

Power consumption: 15W


Weight: 2.5 Lbs

Material: ABS

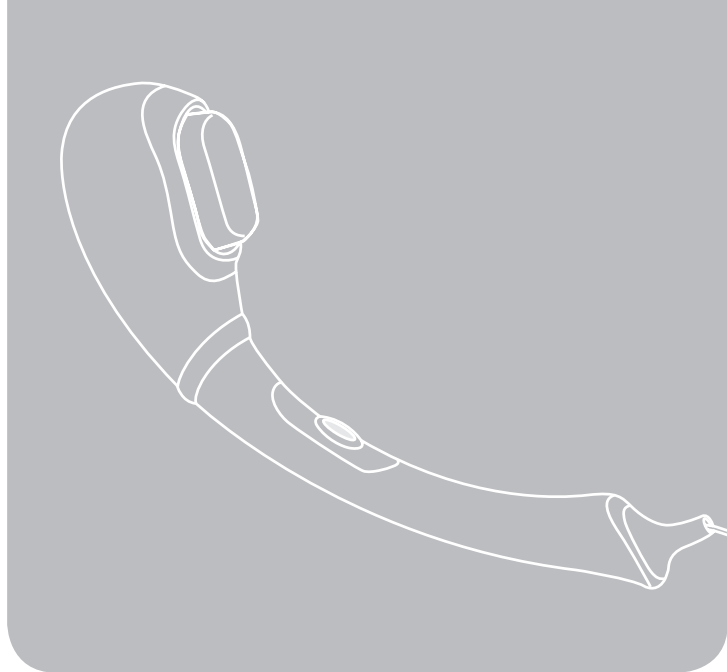**Specifications are subject to change without prior notice.*

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