
humantouch.com
800.355.2762

©2016 Human Touch®, LLC. Patents pending and intellectual property rights enforced.
No medical claims are warranted or implied by the use of this product.
Human Touch is a trademark of Human Touch®, LLC.
HT-ReflexPRO-NA-A0

Register your product **ONLINE** today
www.humantouch.com/register

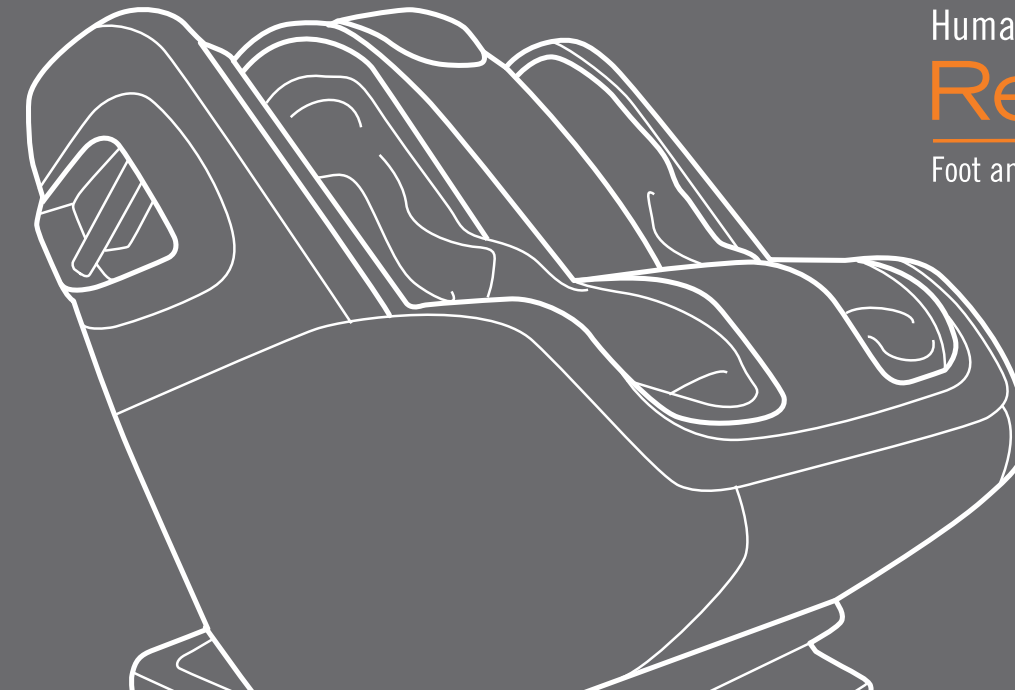Human Touch®
Reflex PRO
Foot and Calf Massager

use & care manual
Human Touch® Reflex PRO

Register your product **ONLINE** today
www.humantouch.com/register

ENG

USE & CARE MANUAL

FRC

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ESN

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For **more than 35 years**, Human Touch has established itself as an industry leader that is changing the lives of its customers with products and services that help them **feel, perform, and live their best**. Developed in cooperation with the medical community and industry experts, we deliver products that blend patented technology, craftsmanship and intelligent ergonomics with superior design to enable therapeutic relief for everyone.

Our industry-leading massage chairs, Perfect Chair® recliners and evolving set of innovative solutions share attributes in their ease-of-use, attention to performance and quality, and design to adapt to most any lifestyle. These attributes are highly valued by our customers and garner broad industry award and recognition. As a result, Human Touch empowers people to enjoy greater health and vitality in their professions, pursuits and everyday lives. And with a population increasingly aware of the benefits of wellness, Human Touch extends its brand leadership and available **solutions for the benefit of everyone**.

Proudly endorsed by:

Register your product at www.humantouch.com to ensure your product is in our system to simplify the warranty process should you ever need it.

humantouch.com

Human Touch® **Reflex PRO** Foot and Calf Massager

welcome to **human touch**®

Your Human Touch® Reflex PRO Foot and Calf Massager is the perfect bridge between state-of-the-art technology and the ancient healing art of massage, emulating the same techniques used by massage professionals. It's a key ingredient for a balanced life — and it's right in your own home!

Our Targeted Relief products deliver therapy to some of the body's most stressed muscles and joints. Beyond delivering unparalleled relaxation and rejuvenation, these amazing devices provide measurable health and wellness benefits that will keep you performing at your best.

The Reflex PRO Foot and Calf Massager incorporates CirQlation®, a life-like massage that delivers the very latest in robotic foot and calf massage technology. We are the only company in the world to use Figure-Eight® Technology, a patented robotic mechanism that massages the calves in an upward, rolling or wave-like motion. This Figure-Eight movement caresses sore calves and rejuvenates tired feet to elevate your relaxation while forcing blood away from the feet and toward the core of the body to improve circulation. Your Reflex PRO incorporates reflexology, a popular holistic and complementary therapy that benefits and improves general health and relieves stress and pain in other parts of the body.

To ensure that you don't overlook any of your foot and calf massager's features, and that you use it safely, read this manual carefully. Then sit back, relax, and enjoy!

keep it safe

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before using this appliance.

 **DANGER: To reduce the risk of electrical shock**

Always unplug this product immediately after using and before cleaning.

 **WARNING: To reduce the risk of burns, fire, electric shock or injury to persons:**

- Do not put hands or fingers near the massage mechanism or rollers while this product is running. Touching the mechanism can result in injury from the rollers squeezing together.
- If the supply cord is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid a hazard
- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Close supervision is necessary when this product is used by, on, or near children, invalids, or disabled persons.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return this product to a service center for examination and repair.
- Do not carry this appliance by supply cord or use cord as a handle.
- Always keep the power cord away from heated surfaces.

- If your feet, calves or legs feel uncomfortable, discontinue use immediately.
- Never operate this product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- Never drop or insert any object into any opening.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- To disconnect, turn all controls to the OFF position, then remove plug from outlet.
- Do not massage any area of the body that is swollen, inflamed or covered with skin eruptions.
- Do not use on unexplained calf pain.
- Do not use this product when standing on a damp floor, or when any part of the body is in contact with plumbing or any similar ground.
- Do not stand on or in appliance. Use only while seated.
- Anyone with a medical condition, including those who are or may be pregnant, should consult a physician before using any massage product.

Safety Precautions

- Comfortable clothing is recommended while using this product.
- Do not sit on the controller.
- Do not pull on the AC cord.
- When using this product after it has been in storage, check to see that the product is operating properly prior to use.

SAVE THESE INSTRUCTIONS

Do not attempt to repair this product yourself.

- This product is designed for personal in-home use. Commercial use will void warranty.
- To reduce the risk of electric shock, this appliance has a polarized plug (one blade is wide than the other). This plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. Do not change the plug in any way.

Care and Maintenance

The product is specially designed and constructed to be maintenance free. No periodic lubrication or servicing is required. Your product should only be serviced by an authorized service center. Foreign objects should not be inserted between the massage nodes or the unit housing. Care should be taken not to overload the unit.

STORAGE:

Your product has been carefully designed to give you many years of trouble-free service. Always unplug the unit when not in use. Additionally, if not in use for a long period of time, we recommend that the power cord be coiled and that the unit be placed in an environment free of dust and moisture. Do not store near heat or open flame. Do not leave this product exposed to direct sunlight for long periods of time. This may cause fading or damage.


 **WARNING:** This appliance is equipped with a thermal resetting protection device. This is a safety feature to protect the appliance from overheating. If the appliance should suddenly stop and will not start, turn off the main switch and pull the plug off the socket outlet. Do not operate the appliance for at least 30 minutes. Failure to turn the appliance off may result in the appliance starting unexpectedly when the device cools.

FABRIC CARE:

To maintain your Human Touch product in optimum condition, dust it frequently and/or vacuum it with a soft brush vacuum extension. Faux Suede products can also be brushed with a suede brush. Wipe with a slightly damp cloth as necessary. For tougher stains on non-leather products, use a mild detergent mixed with water. Leather products should be wiped down occasionally with a leather furniture cleaning solution.

 **WARNING:** Be sure to test any cleaning solution on a small, inconspicuous part of the product to ensure that the color does not fade prior to using the solution on the entire product surface.

get ready

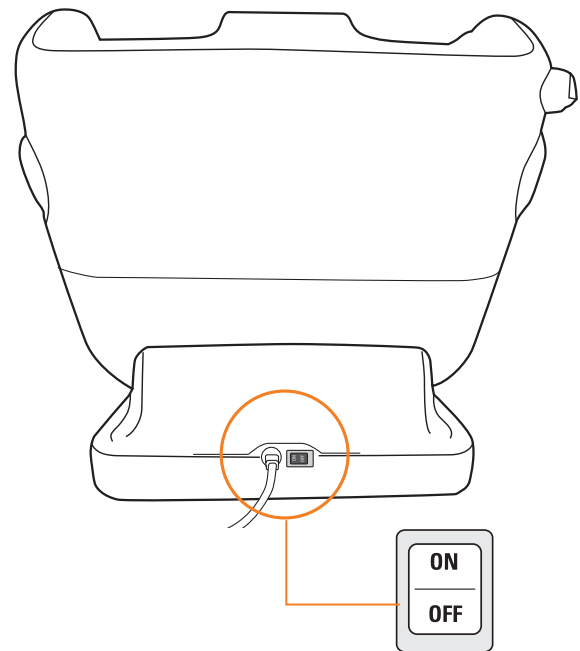1. Plug it in

Plug the power cord into an appropriate power source.

2. Power on


Turn the master power switch (on the rear of the massager) to On.

3. Position your foot and calf massager in front of your favorite chair.

4. Take a seat.

Adjust the massager angle.

Push the angle adjust lever forward, then use your feet to adjust the massager angle as desired. Pull the angle adjust lever backward to lock it into position.

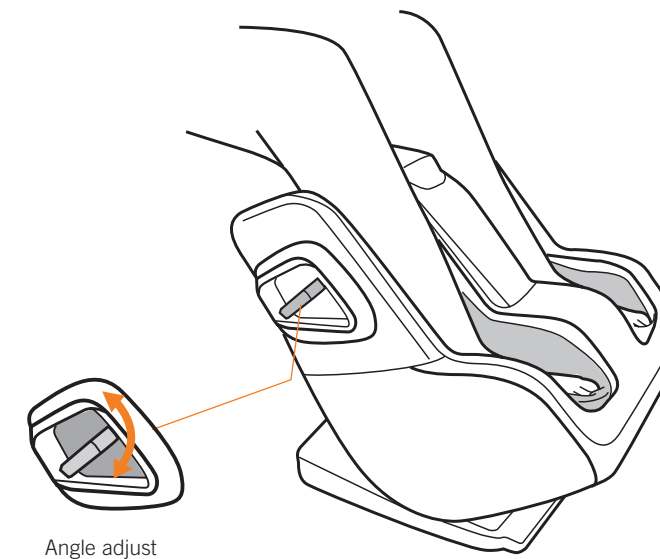Your foot and calf massager comes fully assembled. Simply plug it in, take a seat, adjust the angle, and you're set to go!

NOTE: Follow the below steps to remove your feet from the unit in the event of a power failure:

1. Lift your knees upward (Fig. 1).
2. Rotate your feet sideways and ease them out (Fig. 2).

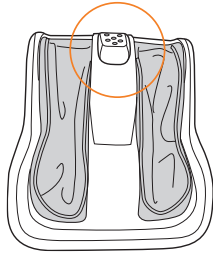
Fig. 1

Fig. 2

 Sitting or standing on your foot and calf massager could result in personal injury and or permanent damage to the massager.

take control

Power – Lights to indicate that the massager is receiving power.

Stop/Park – Returns the massage mechanism to the home position at the bottom of the backrest, then stops.

Vibration – Toggles through two fixed speeds and an auto massage program, which uses varying speeds.

STOP – Cancels the massage and stops all motion.

Manual – Toggles between slow and fast single-speed programs.

Auto – Toggles among four varying-speed programs that incorporate vibration and heat.

Warm Air – Activates foot/calf heat during any massage program.

power stop/park manual relax tone rejuvenate auto refresh auto warm air

human touch®

AUTO PROGRAMS

RELAX – A gentle massage to soothe tired feet and calves.


TONE – A moderate, refreshing foot and calf massage.

REJUVENATE – A firm massage to help sore feet and calves recover.

QUICK-REFRESH – A five-minute medium-to-strong foot and calf massage.


If your feet, calves or legs feel uncomfortable, discontinue use immediately.

Reflexology – Powerful rollers just under the feet provide a robust massage to the sole and heels, providing wellness to your whole body.

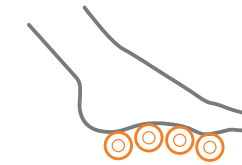Non-Skid Floor Protectors – Floor protectors keep unit in place and ensure quiet operation.

Easy Sleeves®
Removable and washable fabric liners.

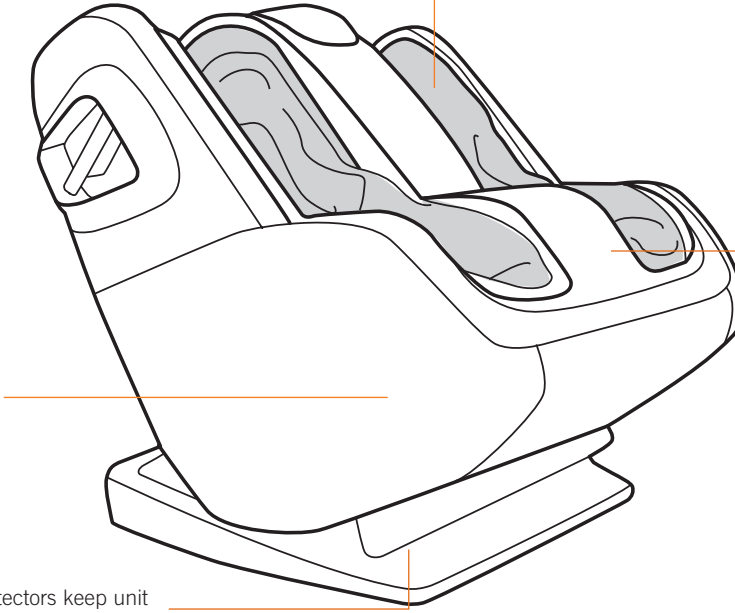Your foot and calf massager is very flexible. You can select from two single-speed or four-multi-speed programs, along with foot/calf vibration and heat!

Patented Figure-Eight® Technology
Refreshes and invigorates sore muscles, while circulating blood from the feet toward the core of the body, delivering the wellness benefits of improved circulation.

Specifications*

Voltage:	AC110-120V 60 Hz
Power Consumption:	0.8 A
Product Weight:	17 lbs.

*Specifications are subject to change without prior notice.

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